**Relaxknow**

# **RXM500**

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Neck and Shoulder Massager

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# DEUTSCH

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# ENGLISH

Technical Data	10
Structure and Accessories	10
Warnings and Precautions	11
Fault Diagnosis	11
Usage	12
Notices / Infos	14
Clean & Maintain	14
Disposal	15

## TECHNICAL DATA

Power: 25 W

Dimensions: L 390 mm × W 340 mm × H 185 mm

Net weight: 2,4 Kg


Voltage: AC100-240 V 50 / 60 Hz - DC12 V, 2.5 A

Scope of delivery: Massage device, AC adapter, car charger, manual

Standard working time: 10 minutes

## STRUCTURE AND ACCESSORIES

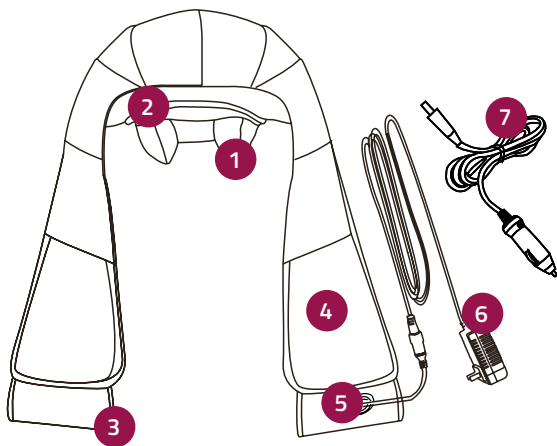**1-** Removable cover (also washable) with zipper. Pleasant net fabric.

**2-** 8 massage balls

**3-** Strap

**4-** PU Leather

**5-** DC Connectors

**6-** Adapter

**7-** Car Charger

## **WARNINGS AND PRECAUTIONS**

Please read these operating instructions carefully before using the device and keep them in a safe place. Use the device only as described and familiarize yourself with all operating and safety instructions before using the product.

- The following people should not use this device: People with osteoporosis, heart pacemakers, heart disease, pregnant women, people with wound injuries and children under 12 years of age.
- **ATTENTION:** In case of acute complaints please consult the doctor before use! The neck massager does not replace any therapeutic massage of a specialist personnel.
- Do not use the unit if it is wet or damp.
- Avoid pulling the power supply unit. If the power supply or adapter is damaged, please stop using it immediately..
- Dismantling or rebuilding the device can cause damage and injury.
- Use the device for approx. 10 minutes and then leave it switched off for approx. 20 minutes.


## **FAULT DIAGNOSIS**

- If the device cannot be switched on, please check that the plug is correctly connected and check the power supply.
- If the noise of the device becomes louder and louder, it may be because it has been in continuous use for too long. Turn off the unit for at least 30 minutes before restarting.
- If the unit turns off during use, it may be because the engine is overheated. Leave the unit switched off for at least 30 minutes before restarting.

**The device is recommended for the following people:**

- People who sit/stand in the same position for a longer time (for example, sitting at a desk in the office for a long time)
- Nervous, tense people
- People with insomnia
- People with stiff neck or stiff shoulder area
- People with headaches or dizziness


You can use the Massage Device for the whole body as shown in the picture:

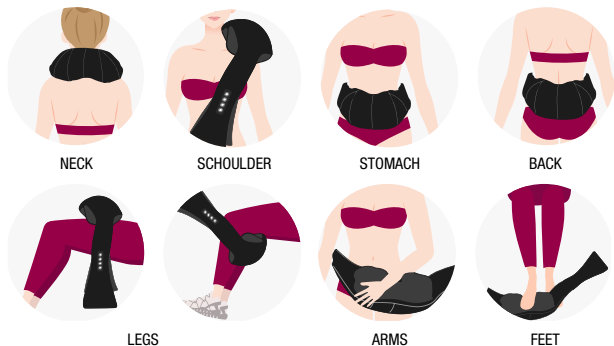Place the massager on the body part that has to be massaged, as shown in the pictures above. The loops, you can put your hands into, are marked with the words left and right. Use these marks to place the massager in the correct position (i.e. keep the left loop on the left and the right loop on the right side). The moving massage elements of the device should be placed on the muscle areas that has to be massaged, i.e. covering the neck, shoulder, abdominal, back, leg, arm, and foot areas. Make sure that the control panel on the left loop is always as close as possible to your hands so that you can manually adjust to the desired setting.

## **USAGE**

Display panel overview:
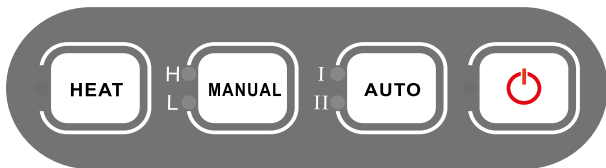
## How to use the Massage Device:

- Plug the AC power cord connector into the main interface of the unit (make sure it is plugged in to prevent it from falling down).
- Plug the AC adapter into the wall outlet. If the indicator light flashes red, the unit is correctly supplied with power.
- Place the Massage Device on the body area of your choice
- Place your hands in the arm straps and keep them taut when applying to the neck, shoulder, back or hip.

**SWITCH ON MESSAGE DEVICE** Press the On/Off button on the control panel. The red control lamp lights up.

**MANUAL MODE** Manual button to adjust the speed and intensity. For a higher speed the „H“ for „Highspeed“ flashes red. For less speed the „L“ for „Lowspeed“ lights up.

**AUTO MODE** Press the „Auto“ button to enter automatic mode. The unit switches off automatically after 10 minutes. Please note that in this mode there is no fluid massage movement, as a natural hand movement imitation takes place.

**HEAT MODE** The massage device is equipped with a warming function. If necessary, proceed as follows: Press the warming function button.

**SWITCH OFF THE MESSAGE DEVICE** Press the On/Off button on the control panel.

## NOTICES / INFOS

- Please wash the removable cover before first usage. Please refer to care instructions in the next section (p. 15)
- The device is for private use only, not for medical or commercial use.
- It is suitable for the whole body massage:
- Neck, shoulder, back, legs, stomach, feet and arms.
- If you have medical concerns, talk to your doctor or health care professional before using the Massage device.
- The massage should always be perceived as pleasant. If you find it uncomfortable or even painful, please change the intensity of the massage and remove some tension from the arm straps.
- Do not use the Massage Device while driving.
- Thanks to the car adapter, you can also use the Massage Device on a long trip as a co-driver or at the rest area, for example.
- After 10 minutes of uninterrupted massage, the device switches off automatically.

## REINIGUNG & PFLEGE

**Attention: Before cleaning the device, make sure that it is switched off and that the AC adapter is disconnected from the wall socket!**

This has to be kept in mind:

- Always use a dry cloth for cleaning.
- Never allow any liquid to penetrate into the unit or accessories as this could cause a short circuit.
- Smaller stains should only be removed with a cloth or damp sponge (with a little liquid mild detergent if necessary).
- Cleaners containing solvents should not be used.
- Do not use the Massage Device until it is completely dry.
- You yourself may only carry out minor cleaning work on the device. Do not repair the unit yourself in case of malfunction. Otherwise, any warranty claim would expire.

Due to the practical zipper, the cover can be easily removed and washed at a temperature of 30°C in the washing machine.

The washable cover is inappropriate for ironing.

The washable cover is inappropriate for bleaching.

## DISPOSAL

Please dispose the package accordingly. Add carton and cardboard to the waste paper and foils in the recyclables collection.

Applicable in the EU and other European countries with systems for separate collection of recyclable materials.

Do not dispose used appliances in household waste. Every consumer is legally obliged to hand in old appliances separately from household waste (e. g. at a collection point). This is the only way to ensure that old appliances are recycled properly and do not have a negative impact on the environment.

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At the following link  
you will find the declaration of conformity:

**[www.relaxnow.de/manuals/pdf/relaxnow\\_rxm500\\_conformity.pdf](http://www.relaxnow.de/manuals/pdf/relaxnow_rxm500_conformity.pdf)**

**Innovamaxx GmbH**

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service@innovamaxx.de

 **+49 30 220 663 551**

[www.relaxxnow.de](http://www.relaxxnow.de)